

Safetec Automatic & Manual Dispensers

Safetec offers an Automatic and Manual Dispenser for Instant Hand Sanitizer; and a Manual Dispenser for SaniWash Antimicrobial Hand Soap. The Automatic Dispenser also has a floor stand option that is easy to place and move around in any facility where hand sanitizer is needed.



**Manual Dispenser
 (Hand Sanitizer & Soap)
 2510037**



**Automatic Dispenser
 (Hand Sanitizer Only)
 27508**



Automatic Dispenser

Drip Tray

**Floor Stand
 (only used with
 Auto Dispenser)**



Front

Back

Drip Tray 2510122

This drip tray can be used on both Automatic and Manual Dispensers

The product also has adhesive placed on the back to easily adhere to surfaces, as an alternative to using screws



Automatic & Manual Dispensers

Instant Hand Sanitizer (Fresh Scent)

- Kills 99.9% of germs (*Active ingredient: 66.5% Ethyl Alcohol*)
- Formula evaporates completely, leaving no sticky residue on hands
- Enriching aloe vera adds moisture with each use, increasing skin hydration

Item#	Description	Case Qty.	Case Wt.
2510037	Manual Dispenser	12 per cs.	tbd
17360	800ml Manual Dispenser Refill	12 per cs.	22 lbs.
27508	Automatic Dispenser	12 per cs.	24 lbs.
17369	800ml Auto Dispenser Refill	12 per cs.	22 lbs.
27507	Floor Stand	1 each	tbd
2510122	Drip Tray	1 each	tbd

Compliance

Meets OSHA, APIC and CDC hand washing recommendations in the absence of soap and water

**Meets ADA operability*



SaniWash® Antibacterial Hand Soap (Manual Dispenser Only)

- Aids in reducing the risk of cross-contamination
- Colorless soap, with a mild pleasant fragrance
- Uses PCMX for antimicrobial action and contains aloe vera for added moisture with each use

Item#	Description	Case Qty.	Case Wt.
2510037	Manual Dispenser	12 per cs.	tbd
34465	800ml Manual Dispenser Refill	12 per cs.	22 lbs
2510122	Drip Tray	1 each	tbd

Compliance

Meets OSHA and APIC hand washing recommendations

**Meets ADA operability*



Industry Uses

Healthcare: Hospitals, physician's offices, labs, pharmacies, long-term care facilities and etc.

Occupational Health & Safety: Government agencies, shipping facilities, schools, jan-san and etc.

Dispenser Instructions

① Choose your Dispenser



**Manual Dispenser
(Hand Sanitizer & Soap)**
2510037



**Automatic Dispenser
(Hand Sanitizer Only)**
27508

② Choose your Refill



**Hand Sanitizer
Manual Refill**
17360



**SaniWash Soap
Manual Refill**
34465



**Hand Sanitizer
Auto Refill**
17369

③ Add a Drip Tray

2510122



Drip Tray

This Drip Tray can be used on both Automatic and Manual Dispensers and is a nice add-on feature.

Adhesive tabs on the back easily adhere to surfaces, as an alternative to using screws.


Manual Dispenser Mounting Instructions

Tape Mounting

1

1. Only recommended for mirrors or other similar flat surfaces
(Clean flat surface with towelette with an abrasive cleanser prior to using towelette. Dry surface thoroughly)
2. Remove backing from tape and press firmly against mounting surface for 10 seconds. Open cover and press directly on back plate for additional 10 seconds. **Allow 8 hours before loading unit.**
3. **Opening Instructions:** Push up on two release buttons on each side of base with thumb and forefinger. With other hand, gently pull cover away from wall at top of dispenser.

(For locking units – insert provided key into keyhole in right-hand release button and press upward; key must remain in keyhole while closing dispenser)

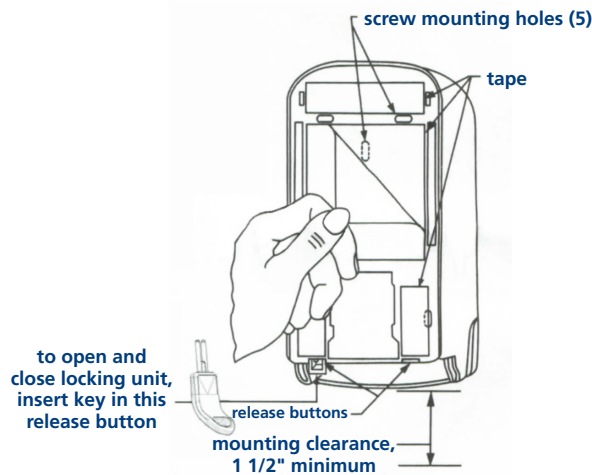

Screw Mounting


1. Recommended for surfaces that are not ideal for Tape Mounting
2. Hold dispenser against wall in desired position, open cover and mark screw hole location on wall, insuring screw heads will not interfere with release bar.
3. Use appropriate anchors and drill size for #10 x 1" self-taping screws to attach dispenser firmly to the wall.

Manual Dispenser Refill Instructions

2

1. Press finger into perforated area and remove key-hole perforated shaped area.
2. Locate nozzle and pull down through bottom of box.
3. Place cartridge on dispenser platform. Pull nozzle with thumb and forefinger to rotate top of tub so Arc of Half Moon Pump is facing toward you.

Arc of Half Moon Pump

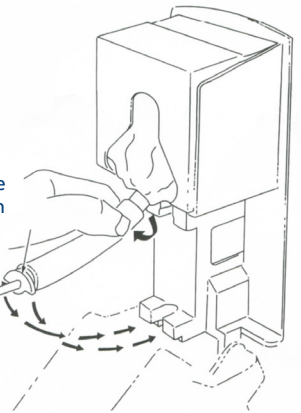


Upper Locator

4. Place top of tube firmly in neck of upper locator and fit Ears of Half Moon Pump into Lower Locator Slots.

Ears of Half Moon Pump

5. Before closing the cover, check to make sure Arc of Half Moon Pump is facing toward you and Ears of Half Moon Pump are firmly seated in slots of lower locator.

Lower Locator Slots

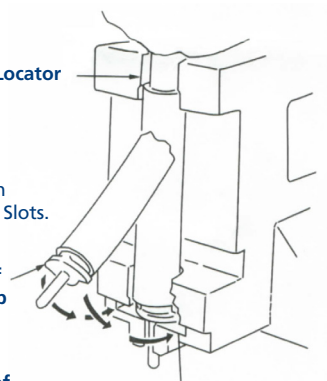


Dispenser Instructions (Continued)

(Automatic Dispenser also available with floor stand option)

Auto Dispenser Mounting Instructions

1 To Open Dispenser

1. Use provided key to press down at top of front cover where key lock is located. Grab the dispenser cover from the top and pull toward you. Cover is hinged at bottom.

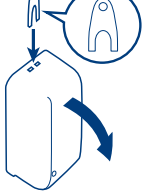

Tape Mounting

1. Only recommended for mirrors or other similar flat surfaces.
2. Dispenser must be mounted at least 10" above any counter or object surfaces.
3. Remove backing from tape and press firmly against mounting surface for 10 seconds. Open cover and press directly on back plate for additional 10 seconds. **Allow 1 hour before loading unit.**

Screw Mounting

1. Hold dispenser against wall in desired location, open cover and mark screw hole location on wall.
2. Install provided anchors into drywall and secure dispenser with provided screws.


Auto Dispenser Set-Up Instructions

2 Installing Batteries (batteries not included)

1. Open dispenser, locate battery compartment on the right side of the dispenser and slide cover forward, away from wall. Place (3) C-cell Alkaline batteries in proper orientation and slide cover back into place.

Installing Hand Sanitizer Refill

1. Press finger into perforated area and remove key-hole perforated shaped area. **Round Disk of Nozzle**
2. Locate nozzle and pull down through bottom of box.
3. Place cartridge on dispenser platform. Pull nozzle with thumb and forefinger, pull round disk of nozzle all the way to bottom of lower locator.
4. Before closing cover, check to make sure round disc of nozzle is firmly in place at bottom of lower locator.
5. After Refill is properly placed, close dispenser cover. LED indicator will blink GREEN/RED 3 times.
6. Place hand underneath dispenser, within 8cm of sensor. LED will blink GREEN and dispense hand sanitizer into hand.




*To Reset: Open and close cover

Sensitivity Adjustment


If sensor is too sensitive, remove black sticker that is covering Sensor Adjustment screw. Turn screw counterclockwise slightly. (i.e. like a dial on a clock, turn in 15 minute intervals to desired setting)


Quick Visual Automatic Installation Instructions

Recommended Installation Height

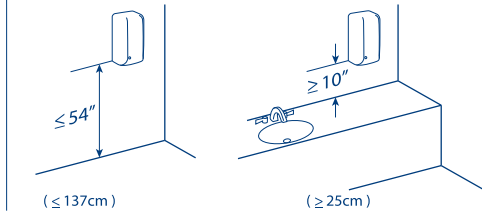

Open Dispenser

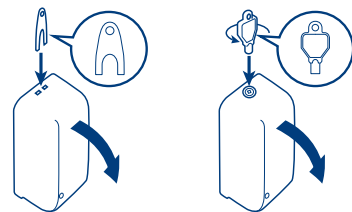


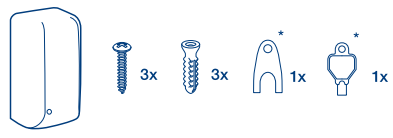

*may not be included in certain models


Installation with Screws

1. Remove the back cover of the dispenser.
2. Insert the screws and anchors into the wall.
3. Attach the dispenser to the wall and secure it with the screws.

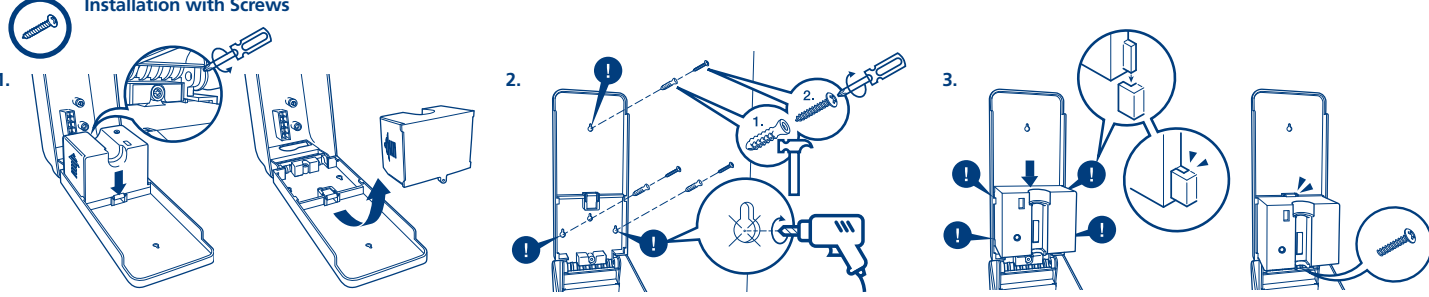


Installation with Adhesive Backing (no screws)

1. Peel off the adhesive backing from the dispenser.
2. Press the dispenser firmly against the wall.

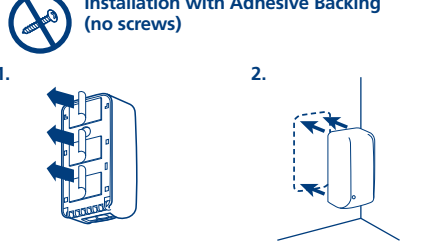


Installation Cartridge

1. Insert the cartridge into the dispenser.
2. Push the cartridge down until it is fully seated.

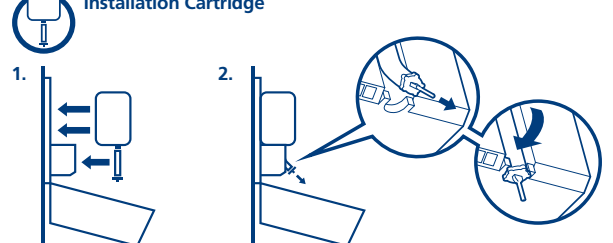


Install Batteries

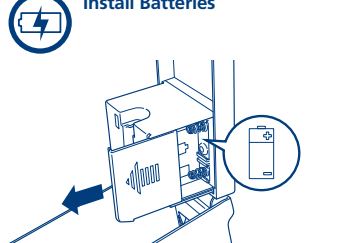


Floor Stand Instructions

Package Include:


- #A. Base Plate (1)
- B. Face Plate (1)
- C. Bottom Pole (1)
- D. Top Pole (1)
- E. Wing nuts (2)
- F. Metal Clips (3)
- G. Plastic Plug (1)

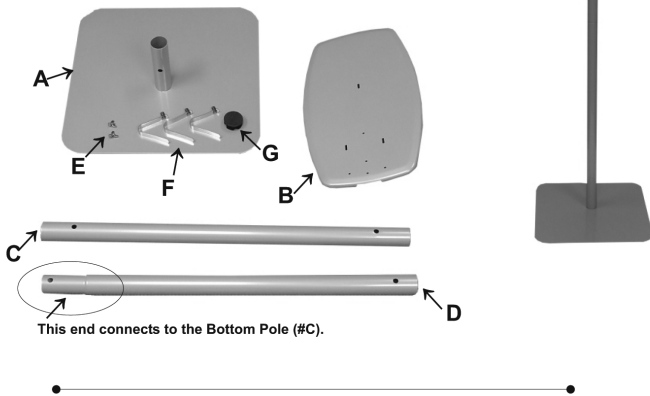


Step 1:
Place a Metal Clip (#F) into the Base Plate (#A).

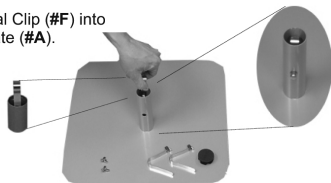


Step 2:
Place a Metal Clip (#F) on each end of the Top Pole (#D).

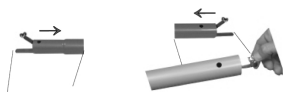


Step 3:
Slide one end of the Bottom Pole (#D) into the Base Plate (#A). Until the **button** pops into place.

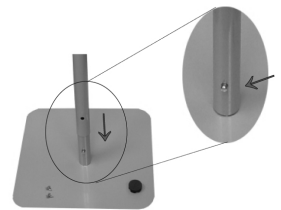

Step 4:
Connect the Top Pole (#D) with the Bottom Pole (#C) by sliding the smaller end of the Top Pole (#D) into the Bottom Pole (#C). Until the button pop into place.

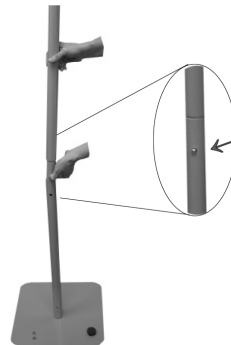

Step 5:
Slide the Face Plate (#B) onto the Top Pole (#D) until the button pops into place.

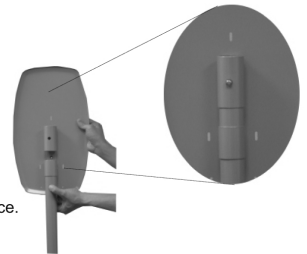

Step 6:
Place the Plastic Plug (#G) onto the Top Pole (#D) by pushing down until flush.

Mounting holes for dispenser

Step 7:
Mount dispenser by using adhesive tape or affix with Wing Nuts (#E) to the Face Plate (#B).

