

MONITORING & DEFIBRILLATION

Product Description	Part Numbers	Page
Atlas™ Monitor	622NP-E1	5
Propaq® CS	01031633	6
Propaq Encore®	01020821	6
Acuity® LT Central Monitoring Station	(various)	7
Welch Allyn AED 10™	970300	8
Welch Allyn AED 20™	972212	8
Welch Allyn PIC 30™	973039COO	9
Welch Allyn PIC 40™	973092	9
Vital Signs Monitor 300 Series	53NTP-E1	10

PRE-OP & RECOVERY

Product Description	Part Numbers	Page
Spot Vital Signs®	42NTB-E1	12
Braun ThermoScan® Pro3000 Ear Thermometer	03000-200	13
SureTemp® Plus 692 Electronic Thermometer	01692-200	13
Universal Desk Charger Set	71641	14
767 Wall Aneroid	7670-01	15
767 Mobile Aneroid with SureTemp® Plus Thermometry	7670-04T9	15
Professional Stethoscope	5079-135	16
Welch Allyn Disposable & Durable Blood Pressure Cuffs	(various)	16
LS-150 Halogen Exam Light	44500	17
CP 20i™ ECG	11103-500	17

SURGERY & PROCEDURES

Product Description	Part Numbers	Page
Fiber Optic Laryngoscope Blades & Handles	(various)	19
Intubating Fiberscope	65301	19
MFI Solarc™ Headlight System	46120	20
MFI Xenon Headlight System	49820	20
SFI / CL300 Headlight System	(various)	21
SFI / CL100 Headlight System	(various)	21
Video Colposcope with Vertical Stand	88000A	22
Exam Light III™	48700	23
KleenSpec® Disposable Vaginal Specula and Illumination System	58001, 78810	23
Electrosurgical System	88530	24
Cryosurgical System	88372	24
Fiber Optic Anoscope/Sigmoidoscope Set	35303	25
LS-200 Procedure Light	44200	25

AMBULATORY SURGERY

4341 State Street Road, P.O. Box 220, Skaneateles Falls, NY 13153-0220 USA
(P) 800-535-6663 (F) 315-685-3361 www.welchallyn.com

WelchAllyn

Advancing Frontline Care™

Welch Allyn thanks St. Joseph's Northeast Surgery Center, Fayetteville, NY, for the generous use of their facilities.
©2003, Welch Allyn, Inc. Printed in USA SM2491 Rev B

Nellcor® is a registered trademark of Nellcor Puritan Bennett Inc. ThermoScan® is a registered trademark of Braun GmbH.

WelchAllyn

Advancing Frontline Care™

Welch Allyn. The Point of Difference at the Point of Care.

It's the point where you and the patient meet. The point of contact, where you have to evaluate, diagnose, monitor, and provide treatment—quickly and accurately, without losing sight of the patient as a person.

The very point, in fact, where we at Welch Allyn focus *our* skills. As the first name in medical diagnostic instruments, we know your ambulatory surgery requirements better than anyone – be it office-based, free standing or hospital-based surgical procedures. The value of time, the need to act quickly, the importance of having absolute confidence in your diagnostic information, the ability to deliver the best care in an increasingly challenging, fast-paced environment.


The result is innovation which puts state-of-the-art technology in the ambulatory surgery setting, with practical, efficient, reliable, affordable instruments that give you the very best tools for an extremely demanding environment.

It's the point where what we do best enhances what *you* do best. The point of care.

AMBULATORY SURGERY

- 1 Atlas™ Monitor #622NP-E1 Pg. 5
- 2 Disposable Blood Pressure Cuffs Pg. 16
- 3 PIC 30™ Monitor/Defibrillator #972039CPO Pg. 9
- 4 Fiber Optic Laryngoscope Handles & Blades Pg. 19
- 5 SureTemp® Plus 692 Electronic Thermometer #01692-200 Pg. 13
- 6 Professional Stethoscope #5079-135 Pg. 16
- 7 Braun ThermoScan® Pro3000 #03000-200 Pg. 13
- 8 Universal Desk Charger Set #71641 Pg. 14
- 9 PanOptic™ Ophthalmoscope Head #11820 Pg. 14
- 10 Micropaq® #01080001 Pg. 7
- 11 Acuity® LT Central Monitoring Station Pg. 7
- 12 Spot Vital Signs® #42NTB-E1 Pg. 12
- 13 Durable One-Piece Blood Pressure Cuffs Pg. 16
- 14 LS-150 Halogen Exam Light #44500 Pg. 17
- 15 CP20i™ ECG Monitor #11103-500 Pg. 17

For more information on the products in this catalog, visit our website at www.welchallyn.com

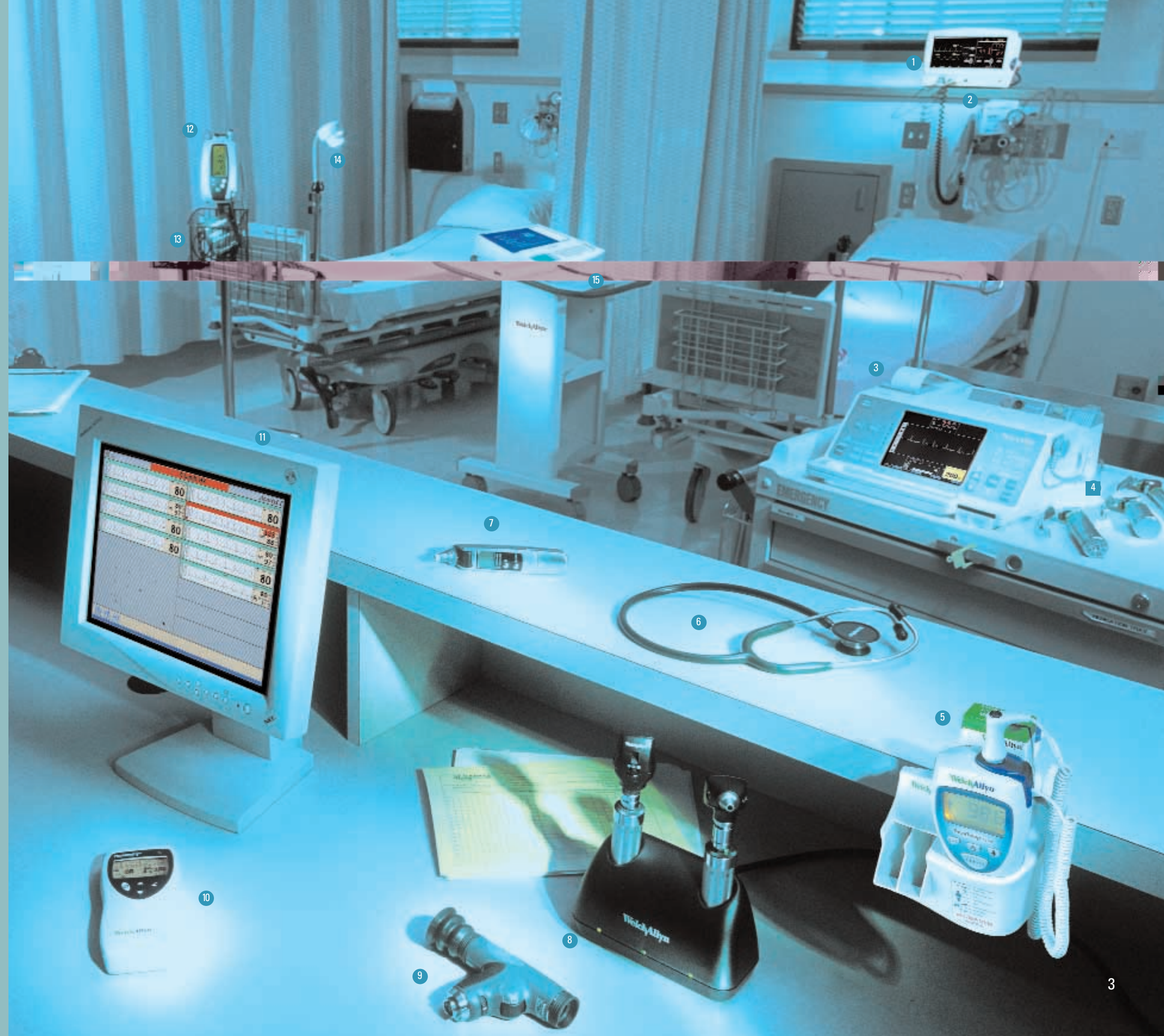

“ I'm impressed with the wide range of cost-effective solutions Welch Allyn offers for conscious sedation and procedure monitoring. Adding a full defibrillation line makes this an even more impressive offering. ”

— TIMOTHY MCCALL, MD-BOARD CERTIFIED ANESTHESIOLOGIST, CAZENOVIA, NEW YORK

need black and white original image


Atlas™ Monitor

A cost-effective monitor with all the essential parameters for procedure rooms & recovery.




- Bright, crisp waveforms and multi-colored LED displays for across-the-room visibility
- Simple, straightforward operation
- Lean, practical, hard-working monitor with configuration options to meet your individual needs
- All monitors include: ECG, SpO₂, and NIBP. Optional parameters include: temperature, respiration, ETCO₂, and printer
- Two-year warranty

- #622NP-E1 Atlas with ECG, Nellcor® SpO₂, NIBP, Temperature, Respiration, Battery, RS232, Printer
- #623NP-E1 Atlas with ECG, Nellcor® SpO₂, NIBP, ETCO₂, Temperature, Respiration, Battery, RS232, and Printer
- #621NP-E1 Atlas with ECG, Nellcor® SpO₂, NIBP, Printer
- #6200-30 Mobile Stand

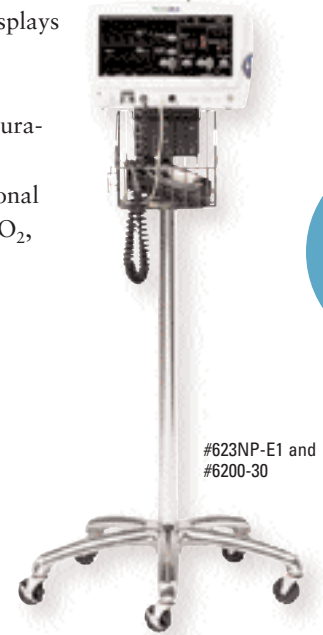

NOW WITH MOTION TOLERANT SpO₂

#623NP-E1

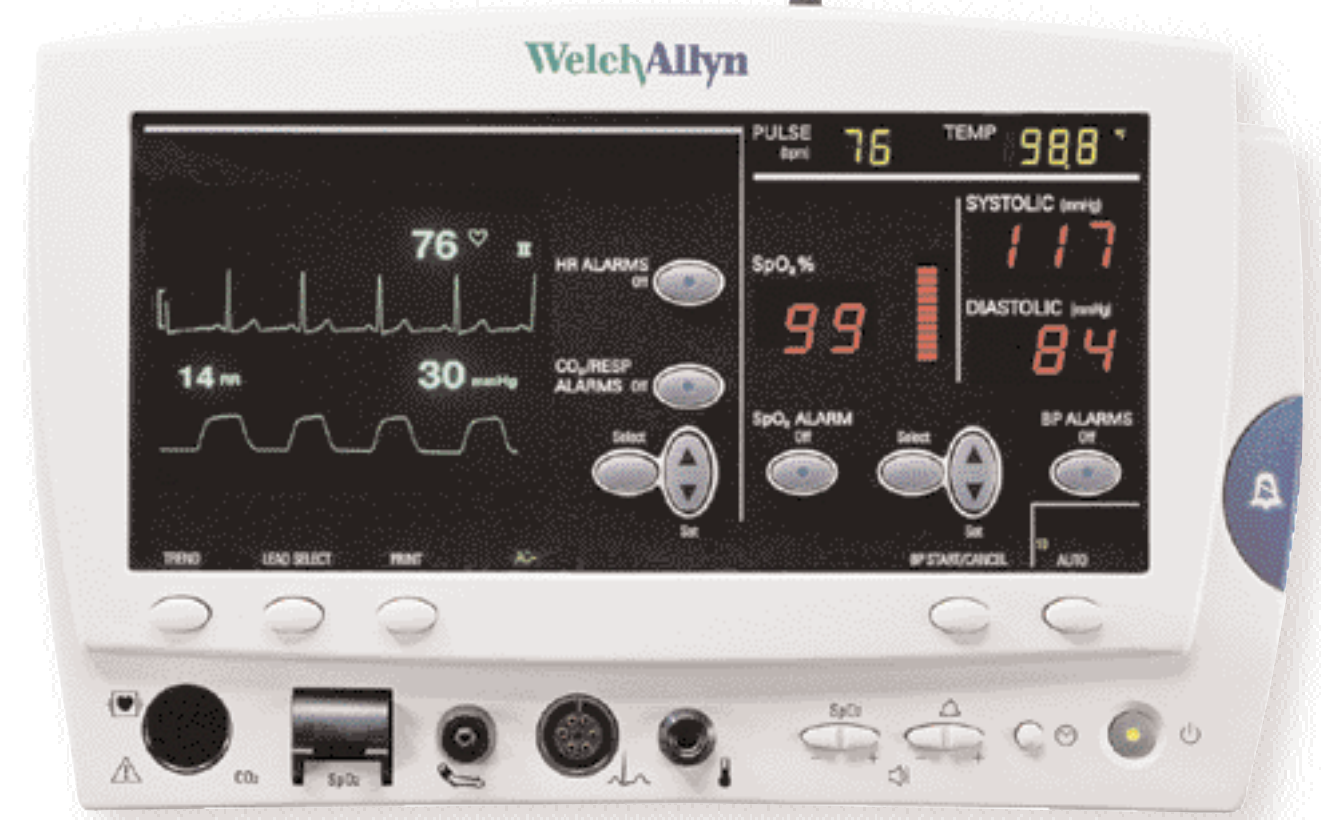


G.I./Endoscopy—
Welch Allyn monitors
automatically down-
load to Provation
automated documenta-
tion, coding and
compliance software.

Propaq® CS Vital Signs Monitor

Small, light, and easy to use, it's the perfect instrument for advanced vital sign monitoring in all ambulatory surgery environments.

- Large color display may be read across the room with color-differentiated parameters
- Custom selection of up to four viewable waveforms
- Adult, pediatric, and neonatal modes for customized monitoring of specific patient populations
- Optional—Wireless connectivity to the Acuity and Acuity LT Central Station allows truly flexible monitoring
- Three-year warranty

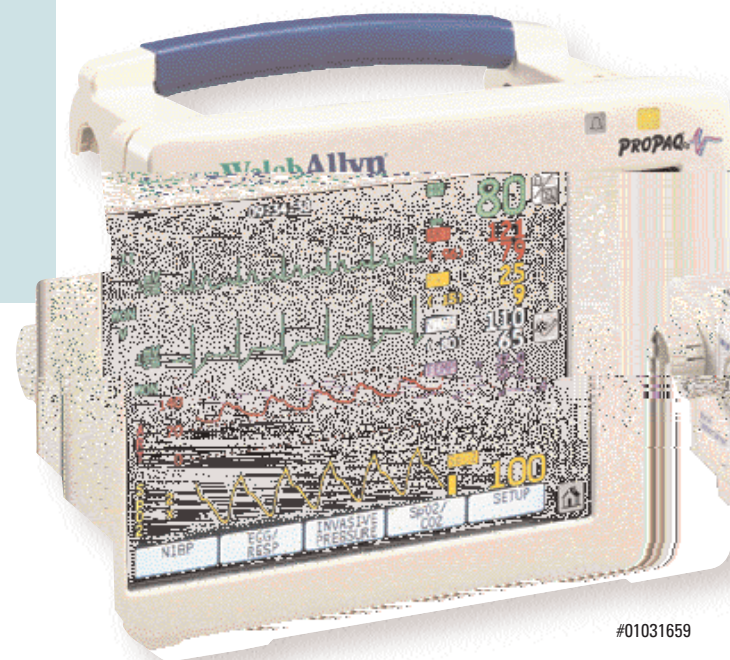

#01031659

- #01031633 Propaq CS, Model 242, Monitors 3 Lead ECG, Respiration, NIBP, Motion-Tolerant SpO₂, and 2 Channels of Temperature
- #01031659 Propaq CS, Model 246, Monitors 3 Lead ECG, Respiration, NIBP, 2 Invasive BP, Motion-Tolerant SpO₂, Expansion Module with Printer, 2 Channels of Temperature, and Dualstream CO₂

Acuity® LT Central Monitoring Station

The affordable, self-deployable central monitoring solution that receives, displays, and stores vital signs information for up to 12 patients from a wireless or hardwired system of monitoring devices.

- Receives, displays and stores vital sign information for up to 12 patients from either wireless or hardwired systems of patient monitors
- Remote control of ambulatory or bedside patient monitors from a central station
- Twenty-four hour full disclosure and life-threatening arrhythmia analysis (Asystole, VFib and VTach alarms)
- Simplified installation allows your facility the cost-saving option of completing your own installation
- Self-paced clinical training CD included
- Use in conjunction with Welch Allyn's Propaq and Micropaq ambulatory monitors




Contact your local Welch Allyn sales representative to plan the ideal Acuity LT configuration for your needs

Propaq Encore® Monitor

Rugged enough for the fast-paced monitoring demands of Ambulatory departments and patient transport

- The gold standard for transport monitoring
- Long battery life protects against data loss during transport and procedures
- Lightweight and durable
- Motion-tolerant SpO₂ and NIBP for accurate, consistent readings in challenging environments
- Low cost of ownership

- #01020821 Propaq Encore, Model 202-EL, Monitors 3-Lead ECG, Respiration, NIBP, Motion-Tolerant SpO₂, and 2 Channels of Temperature
- #01020846 Propaq Encore, Model 206-EL, Monitors 3-Lead ECG, Respiration, NIBP, 2 Invasive BP, Motion-Tolerant SpO₂, Expansion Module with Printer, 2 Channels of Temperature, and Dualstream CO₂

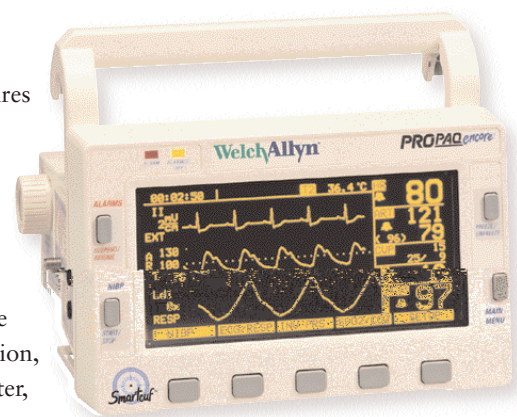

#01020846

#01080005



Micropaq®

The convenience and efficiency of a patient monitor that's wireless and wearable. Ideal for ambulating patients and administering medications at the bedside

- Displays heart rate, one or two ECG leads, motion-tolerant SpO₂ numerics and pulse bar
- Backlit LCD with automatic power save mode for night time use
- Continues monitoring even when out of telemetry range
- Up to 25-hour battery life when fully charged
- Connects wirelessly to Welch Allyn's Acuity and Acuity LT Central Monitoring Stations
- Alarms may be suspended and silenced locally at the monitor, or remotely at the central station

- #01080004 Micropaq, ECG only
- #01080005 Micropaq, ECG and Masimo® Pulse Oximetry
- #01080006 Micropaq, ECG and Nellcor® Pulse Oximetry





#972214

Welch Allyn AED 20™

The Welch Allyn AED 20 provides basic AED functionality with optional Lead II monitoring using snap-based electrodes. Optional manual defibrillation (2j to 360j) makes this the first choice for facilities with multiple levels of first responders.

- Choice of rechargeable battery for frequent users or a long life non-rechargeable battery for infrequent rescuers
- Escalating energy (200j, 300j, 360j), per AHA guidelines
- Energy protocol can be changed to meet specific needs
- Optional ECG monitoring and optional passcode protected manual defibrillation
- Data card transfer and direct basic print capability of critical patient data

- #972212 Welch Allyn AED 20 with Rechargeable NiMH Battery and Charging Station
- #972214 Welch Allyn AED 20 (Rechargeable) with ECG Waveform and Single Channel Monitoring
- #972216 Welch Allyn AED 20 (Rechargeable) with ECG, Monitoring, and Manual Defibrillation

Welch Allyn AED 10™

Designed for easy operation and quick deployment by the infrequent user, the AED 10 is the solution for quick and easy response by minimally trained personnel.

- Text and voice instructions guide the user through the rescue
- Color-coded pads and illustrations ensure proper pad placement
- Escalating energy (200j, 300j, 360j), per AHA guidelines
- Energy protocol can be changed to meet specific needs
- High performance battery delivers 5 years of standby performance or over 90 discharges at 360 joules
- Internal memory stores event data including ECG waveforms, time / date stamps, shock delivery and system actions
- Compatible with major manufacturer defibrillator pads

- #970300 Welch Allyn AED 10
- #970302 Welch Allyn AED 10 with Soft Shell Carry Case



#970300

Welch Allyn PIC 40™

A multi-parameter monitor, AED, and defibrillator with external pacing, the PIC 40 is the perfect solution for monitoring/defibrillation in ambulatory surgery.

- Radiant color display provides waveforms, patient data and defibrillation functions that can be viewed from virtually any angle
- 3 or 5 lead ECG and SpO₂
- Optional parameters include NIBP and external pacing
- Manual defibrillator provides full range of energy settings (2j to 360j)
- AED is configurable and follows AHA escalating energy (200j, 300j, 360j) guidelines

- #973092 PIC 40
- #973093 PIC 40 with NIBP
- #973094 PIC 40 with NIBP and external pacing
- #973095 PIC 40 with external Pacing




#973094

Welch Allyn PIC 30™

Designed for the crash-cart and other manual defibrillation environments, the PIC 30 monitor/defibrillator is the economical choice for hospitals and surgery centers.

- Radiant color display provides excellent visibility
- Annotating recorder captures all vital data on 50mm printer
- Full manual defibrillation (2j to 360j) and sync capabilities
- Optional external pacing and shock advisory
- Also available with monochrome display

- #973039C00 PIC 30 Color Screen
- #973039CP0 PIC 30 with External Pacing
- #973039COA PIC 30 with Shock Advisory
- #973039CPA PIC 30 with External Pacing and Advisory

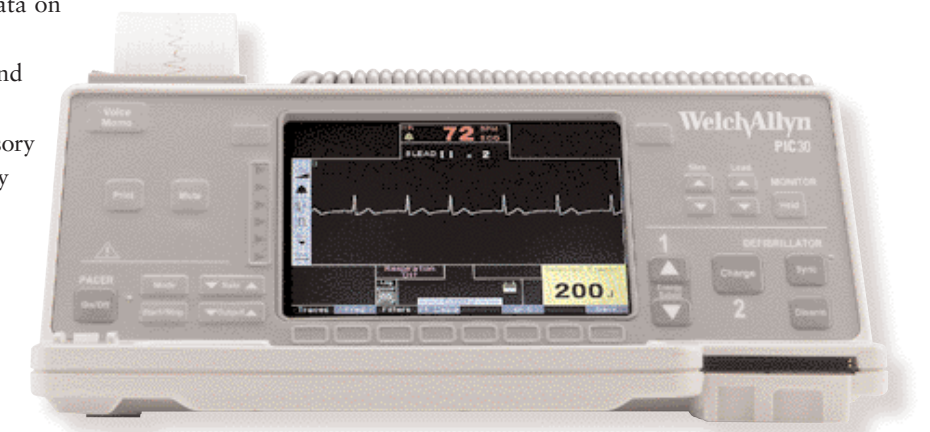

#973039C00

Ask your local Welch Allyn representative how Welch Allyn Vital Signs Monitor and Spot Vital Signs® devices can interface with Electronic Medical Records (EMRs).

Welch Allyn® Vital Signs Monitor 300 Series

For simple and accurate spot-check vital signs assessment or continuous monitoring, the newly redesigned Vital Signs Monitor is the definition of “advanced monitoring made affordable.”

- Standard non-invasive blood pressure (NIBP), pulse rate and MAP
- Optional – SureTemp® Thermometry*, motion-tolerant Nellcor® Pulse Oximetry and integrated thermal printer
- Supports adult, pediatric and neonatal patients
- Simple to use – new LCD display facilitates menu navigation
- Includes automatic blood pressure mode, programmable alarms, built-in memory capability, and intuitive icons
- Lightweight and portable, with a variety of mounting options

#53NTP-E1 Vital Signs Monitor complete with NIBP, Nellcor® Pulse Oximetry, SureTemp® Thermometry and Printer

#53000-E1 Vital Signs Monitor with NIBP (includes Pulse Rate and MAP)

#530T0-E1 Vital Signs Monitor with NIBP and SureTemp® Thermometry

#53N00-E1 Vital Signs Monitor with NIBP and Nellcor® Pulse Oximetry

#53NT0-E1 Vital Signs Monitor with NIBP and Nellcor® Pulse Oximetry, and SureTemp® Thermometry

#008-0835-00 Mobile Stand

*See page 13 for full details on Welch Allyn SureTemp Thermometry

NOW WITH MOTION TOLERANT SpO₂



PRE-OP AND RECOVERY

Vital signs collection is repeated so frequently that it can have a huge impact on your efficiency. Welch Allyn offers many fast, dependable solutions to simplify this process.

— ROBERT CORONA, D.O., F.A.C.P., CHIEF MEDICAL OFFICER, WELCH ALLYN, INC.

“We love the convenience. Being able to automatically collect all the vital signs in a short amount of time really frees up the nurses.”

— LAURIE MCGRAW, RN,
MINONG, WI



Spot Vital Signs®

Take vital signs in seconds with superior technology that spot checks easily, accurately and affordably.

- Automated measurement of non-invasive blood pressure (NIBP), pulse rate, temperature, and Pulse Oximetry in 30 seconds
- All configurations include: NIBP, Pulse Rate, and MAP
- Optional—SureTemp® Thermometry or Nellcor® Pulse Oximetry
- Flexible configurations yet fully upgradeable to meet your future needs
- Lightweight and portable with a variety of mounting options

#42NTB-E1 Spot Vital Signs with NIBP, Nellcor® Pulse Oximetry, SureTemp® Thermometry

#4200B-E1 Spot Vital Signs with NIBP

#420TB-E1 Spot Vital Signs with NIBP and SureTemp® Thermometry

#42NOB-E1 Spot Vital Signs with NIBP, Nellcor® Pulse Oximetry

#4200-60 Mobile Stand



#42NTB-E1 and
#4200-60

To meet the needs of a busy environment, Welch Allyn offers you a variety of modalities for temperature-taking for all types of patients.

BRAUN ThermoScan. Pro3000 Electronic Ear Thermometer


For patients and situations that require fast action, take temperatures in a second—literally—with infrared technology that measures body temperature accurately through the ear.

- Sensor scans eight temperatures in one second, displaying the highest reading
- Built-in safeguards virtually eliminate risk of cross-contamination
- Soft, comfortable probe cover keeps probe tip clean for accuracy
- Easy-to-read LCD in Fahrenheit or Celsius
- Takes approximately 5,000 readings on two lithium batteries

#03000-200 Braun ThermoScan® Pro3000


#03000-300 Braun ThermoScan® Pro3000 with Wall Mount

#05074-005 Probe Covers (Box of 5000)

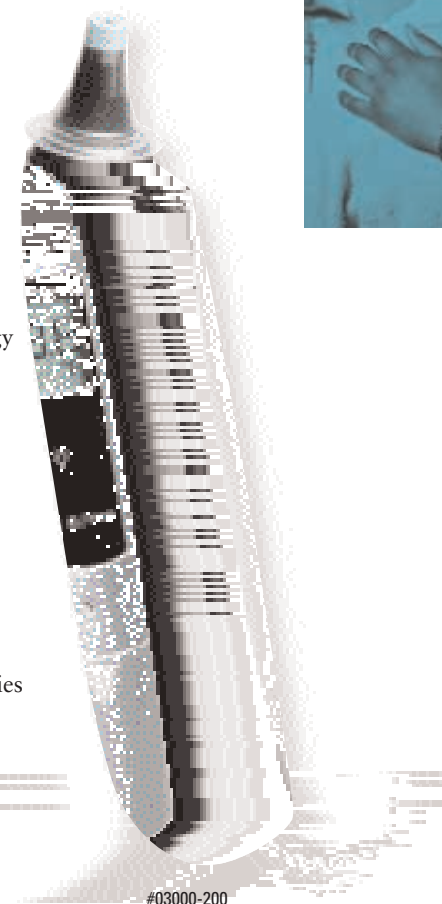

#03000-200

NOW
WITH MOTION
TOLERANT
SpO₂

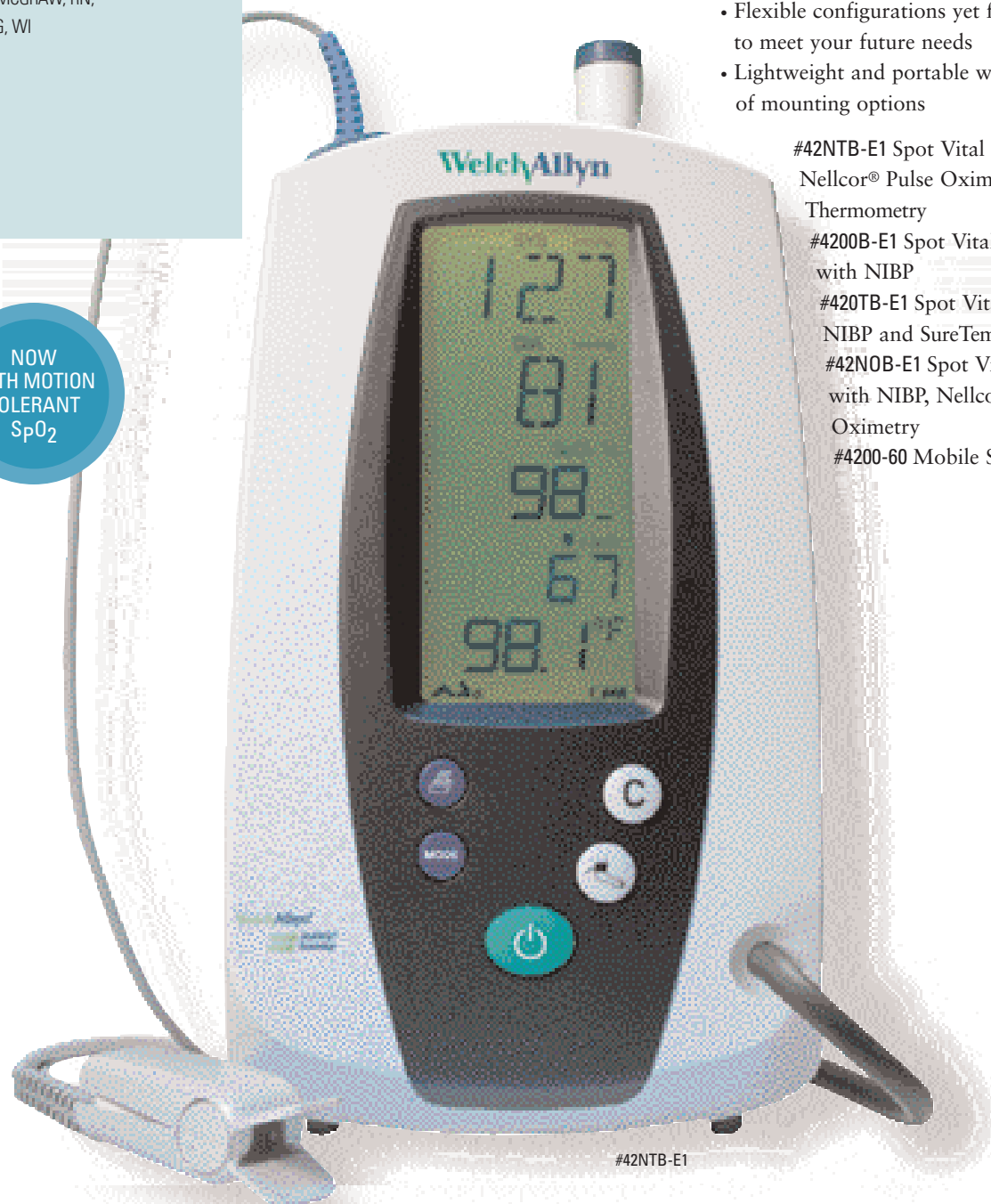

#42NTB-E1

SureTemp® Plus 692 Thermometer

Welch Allyn's newly redesigned SureTemp Plus combines the speed of ear thermometry with the accuracy and reliability of traditional oral thermometry. You can take fast, four-second oral, 10-second pediatric axillary, 15-second adult axillary and 10-second rectal temperatures.

- New ergonomic design – comfortable to hold, easy to use, and very durable
- Interchangeable, removable probe wells to reduce the risk of cross-contamination
- Large backlit LCD display with location identification capability
- Includes a 60-second pulse timer for counting pulses and respirations
- Your choice of security alarm systems to help prevent theft
- Takes approximately 5000 readings on 3 AA batteries


#01692-200 SureTemp Plus 692 Thermometer, Oral Probe with Well, Wall Mount and Security System

#01692-700 SureTemp Plus 692 Thermometer, Oral Probe with Well, Wall Mount and Rolling Stand

#02892-000 Rectal Probe and Well Kit (4 ft)

#02893-000 Oral Probe and Well Kit (4 ft)

#05031-110 Thermometer Probe Covers (box of 10,000)

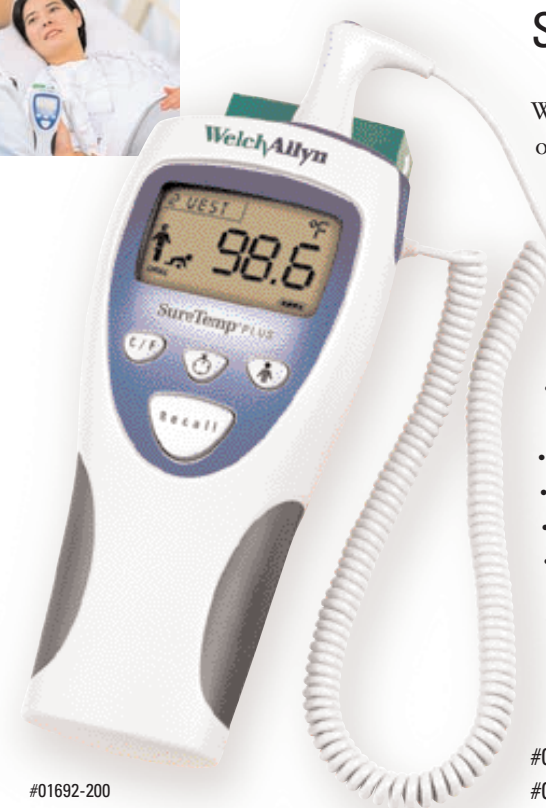

#01692-200

Universal Desk Charger Set

Keep your diagnostic instruments powered up and ready to go with our most popular handles and instruments in one easy-to-order set.

- Colored LEDs signal power is “on” and instruments are charging
- Advanced circuitry ensures constant, safe charging and optimum, long-lasting lamp life
- Optional—wall mount kit and/or specula tray


#71811-PSM
 Universal Desk Charger plus:
 PanOptic™ Ophthalmoscope (#11810)
 Diagnostic Otoscope (#25020)
 (2) Lithium Ion Handles (#71900)

#71811-P
 Universal Desk Charger plus:
 PanOptic™ Ophthalmoscope (#11810)
 Diagnostic Otoscope (#25020)
 (2) Ni-Cad Handles (#71670)

#71641
 Universal Desk Charger plus:
 Coaxial Ophthalmoscope (#11720)
 Diagnostic Otoscope (#25020)
 (2) Ni-Cad Handles (#71670)



#71811-PSM

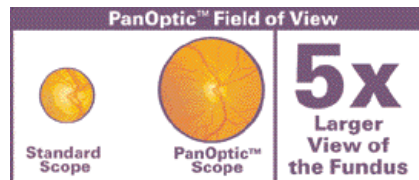

#11720

3.5v Coaxial Ophthalmoscope

- Halogen HPX™ lamp provides 30% whiter, brighter light for true tissue color and long performance
- Detect corneal abrasions and foreign bodies with cobalt blue filter
- Coaxial optical system produces shadow-free spot viewing
- Polarizing filter virtually eliminates corneal reflection
- Sealed optics keep out dust and dirt

#11720 3.5v Coaxial Ophthalmoscope (Head Only)


Tycos® 767-Series Wall and Mobile Aneroids

The convenience of blood pressure at every bed. Clinical-grade accuracy and mercury-free design make the 767-Series one of the best Tycos® aneroids ever for manual blood pressure.

- Certified accurate to ± 3 mmHg
- Mercury-free and latex-free for safety
- Ruby jeweled movement for long life
- Laser-engraved dial face for accuracy
- Large, easy-to-read dial face (6.25" / 15 cm)
- Convenient built-in cuff storage, adjustable-height, 360° swivel stand for clear viewing
- Lifetime warranty

#7670-01 Wall Aneroid with Durable One-Piece Cuff (Adult)

#7670-04 with Four-Leg Mobile Stand, Durable One-Piece Cuff (Adult), and Base Weight

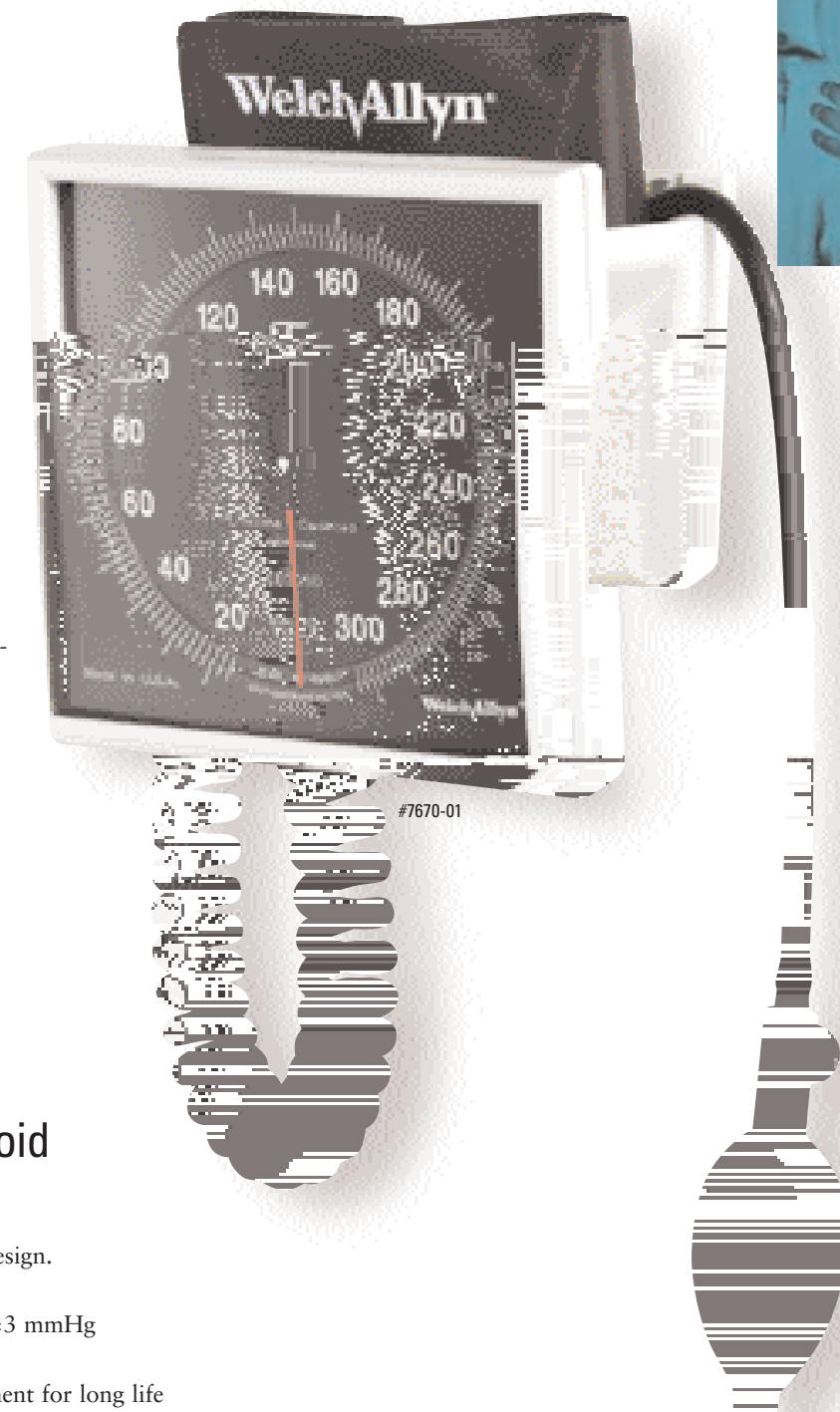

#7670-01


Tycos® Classic Hand Aneroid

Unsurpassed reliability in a traditional, handheld design.

- Certified accuracy to ± 3 mmHg
- Latex-free for safety
- Ruby jeweled movement for long life
- Laser-engraved dial face for accuracy
- Durable one-piece cuff
- Lifetime warranty

#5098-02 Hand Aneroid with Durable One-Piece Cuff (Adult)

#5098-28 Hand Aneroid with Multi-cuff Kit (Adult, Large Adult and Child Print Cuff)

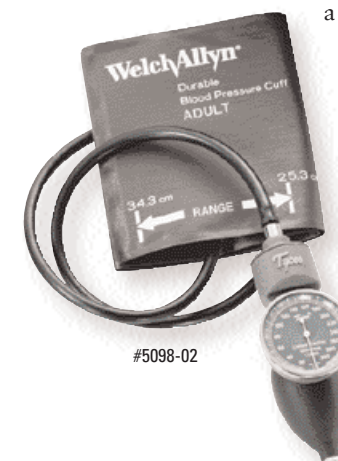

#5098-02

WELCH ALLYN OFFERS FREE MERCURY DISPOSAL. Mercury spills are a costly and hazardous distraction. In fact, some parts of the country are regulating and even banning medical products containing mercury, such as thermometers and sphygmomanometers. Welch Alllyn is helping—with our own mercury-free sphygmomanometers, and with a free program that will safely dispose of your existing mercury sphygmomanometers as part of a facility-wide conversion program. For more information, contact your local Welch Alllyn representative or visit the EPA website at <http://www.epa.gov/seahome/mercury/src/outmerc.htm>



Professional Stethoscope

Excellent acoustics and superior design make this full-range instrument the first choice for professionals who need auscultation accuracy and all-day comfort.

- Meets all diagnostic needs, from detecting Korotkoff sounds of blood pressure to analyzing subtle heart murmurs
- Stainless steel chestpiece intensifies body sounds transmitted through tissues

#5079-135 Professional Adult, Double-head Stethoscope, 28" Tubing
Shown in black. Available in a variety of colors.



Disposable Uniscope

The best protection against cross-contamination.

- Reduces risk of cross-contamination with the accuracy of conventional stethoscopes
- The safest choice for isolation, infectious disease and heavy fluid cases

#17462 Adult, Single-Head Disposable Stethoscope, Dark Blue
#17462P Pediatric, Single-Head Disposable Stethoscope, Light Blue
Available in a variety of colors and sizes.




Lightweight Stethoscopes

Economical and comfortable for general use.

- Quality stethoscope designed for blood pressure measurement and general auscultation
- Lightweight and flexible single-channel tubing provides comfortable use over extended periods

#5079-73 Lightweight Double-Head Stethoscope, Velvet Black
Available in a variety of colors.

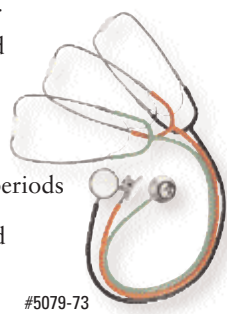

Disposable One-Piece Cuffs and Inflation Systems

Prevent cross-contamination in high-risk areas.

Durable One-Piece Cuffs and Inflation Systems

Never remove, clean or replace a bladder.

- Cuffs are light, durable and latex-free
- Available in one- and two-tube designs for electronic NIBP monitors and sphygmomanometers
- Connectors attached to fit virtually every brand of NIBP monitors

*See page 27 for complete listing of cuff sizes and part numbers

#44500

LS-150 Examination Light

Fifty watts of bright, color-corrected light for easier, faster exams.

- Color-corrected 50-watt lamp delivers significantly whiter light for better tissue diagnosis
- Luminaire stays cool to the touch, even during multiple exams
- Design allows clear view and one-hand positioning
- Caster base allows easy mobility

#44500 LS-150 with Floor Stand
#44510 LS-150 with Wall Mount
#44350 Caster Base

LS-135 Examination Light

High-quality, low-cost halogen illumination.

- Long-lasting, 35-watt halogen lamp produces whiter, brighter light than conventional gooseneck units
- Dual-intensity control switch allows you to instantly adjust light output to match your exact requirements
- Compact base frees up valuable floor space


#44300 LS-135 with Floor Stand
#44310 LS-135 with Wall Mount
#44350 Caster Base

Welch Allyn CP 20™

Part of Welch Allyn's new family of multi-channel electrocardiographs, the full-featured CP 20 is a 12 lead resting ECG that features a full alphanumeric keyboard for fast and easy entry of patient demographics, user-programmable report formats printed on regular, 8 1/2" x 11" thermal paper, ECG measurements, advanced filters for optimal ECG trace quality and battery or AC operation.

- Integrated display to preview the ECG waveform
- Storage for up to 50 ECGs
- Optional interpretation program (adult and pediatric)
- Optional communication option for connectivity to the CardioPerfect™
- Workstation for editing, review and paperless storage
- Three year warranty

#11103-500 Welch Allyn CP 20i ECG with Interpretation
#11103-510 Welch Allyn CP 20i ECG with Interpretation, with Communication Only
#11103-200 Welch Allyn CP 20 ECG without Interpretation
#11103-210 Welch Allyn CP 20 ECG without Interpretation, with Communication Only
#11103-200 Welch Allyn CP 20 Cart with Cable Arm

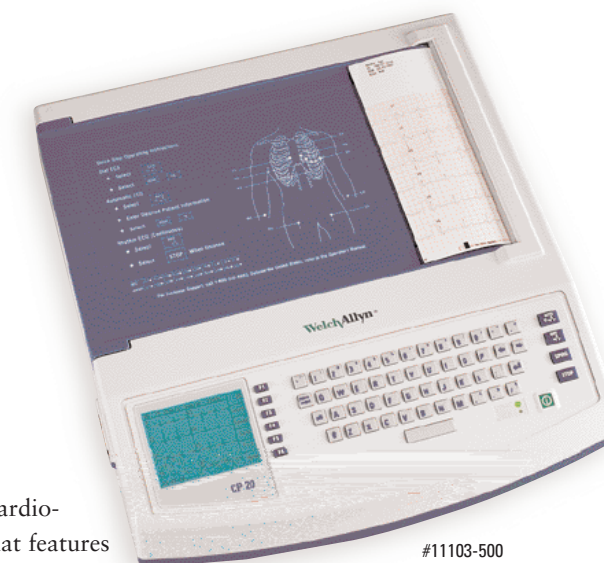

#11103-500

“Modern technology is introducing simpler, less expensive methods of diagnosing and treating various disease entities, which leads to a tremendous savings for health care.

Welch Allyn continues to be in the forefront of those who are developing innovative new equipment to support these efforts.”

— JOHN L. PFENNINGER, MD, PRESIDENT, NATIONAL PROCEDURES INSTITUTE
 CLINICAL PROFESSOR OF FAMILY MEDICINE, MICHIGAN STATE UNIVERSITY
 EDITOR, PROCEDURES FOR PRIMARY CARE



Fiber Optic Laryngoscope Handles

Bright, white HPX™ illumination offers a clear, color-corrected view.

- Patented HPX™ technology blends Halogen and Xenon gases under pressure to produce bright, white light for true tissue color and consistent, long-lasting illumination
- Single lamp cartridge can be used across all sizes and handles
- Up to 40% lighter than our previous design, for increased responsiveness and maneuverability
- Knurled finish ensures secure grip

LIGHTWEIGHT BATTERY-POWERED FIBER OPTIC HANDLES

- #60813 Medium, uses two C-size batteries
- #60814 Penlight, uses two AA-size batteries
- #60815 Stubby, uses two AA-size batteries

Fiber Optic Laryngoscope Blades

- Removable, repolishable light pipes and free repolishing service help blades last longer
- Damaged light pathways can be repaired at a fraction of the cost of replacement blade
- Single-piece stainless steel construction ensures blade integrity

MACINTOSH FIBER OPTIC BLADES

Model	Blade Size
#69061	#1
#69062	#2
#69063	#3
#69064	#4

ENGLISH MACINTOSH FIBER OPTIC BLADES

Model	Blade Size
#69211	#1
#69212	#2
#69213	#3
#69214	#4

MILLER FIBER OPTIC BLADES


Model	Blade Size
#68065	#00
#68060	#0
#68061	#1
#68062	#2
#68063	#3
#68064	#4

Intubating Fiberscope

The ideal instrument for difficult intubations.

- Patented Solarc™ light source delivers whiter, brighter light in a smaller, space-saving package
- Designed for easier, more efficient endotracheal intubation
- Comes with 60 cm flexible probe, 35 mm insertion tube, 1.4 mm multipurpose channel, 90° field of view and 3-50 mm focal range, 260° tip deflection and 25 mm angulation section, and +2 to -8 diopter adjustment

#65301 FL-100 Intubating Fiberscope with Solarc™ Light Source

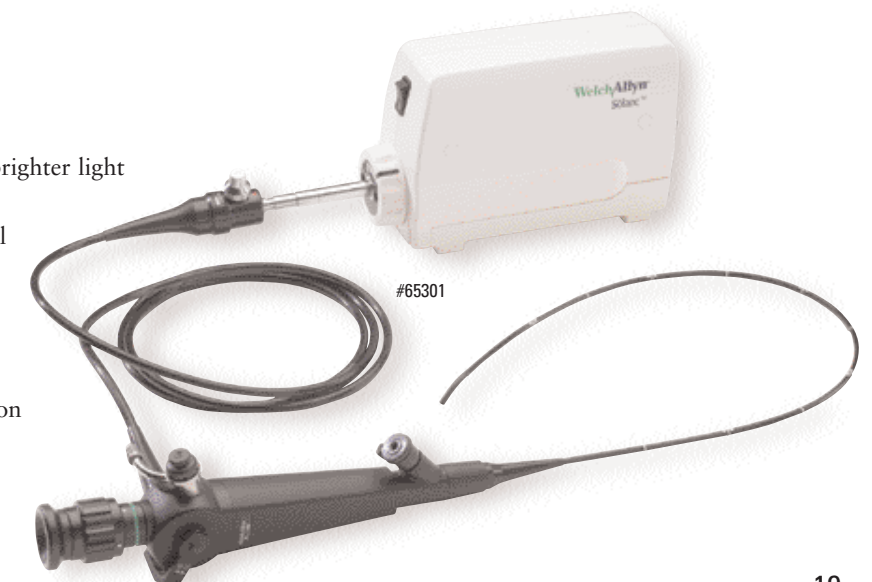

MFI Headlight Systems

Shielded fiber optic light guide and milled aluminum headlight results in unsurpassed durability and stands up to years of demanding use. Greatly reduces ongoing expense of fiber optic bundle replacement common with all other traditional headlight systems.



MFI Solarc™ Headlight System

Xenon quality light at a halogen price

- Available multi-connection turret allows attachment of most Storz, Olympus, Wolf and ACMI fiber optic cables
- Low cost Solarc™ lamp provides on average over 550 hours of illumination
- Compact Solarc™ light source includes mounting bracket for attachment to any flat horizontal surface

- #46120 MFI Solarc™ Headlight System
- #46120TW MFI Solarc™ Headlight System with Turret
- #46100 Solarc™ Light Source Only (ACMI port)
- #46100TW Solarc™ Light Source Only with Turret (ACMI, Storz, Olympus, Wolf)
- #46150 Solarc™ Mobile Stand

MFI Xenon Headlight System

Intense illumination for the most demanding procedures.

- 300-Watt Xenon light source provides intense illumination
- Xenon lamp module provides on average over 500 hours of illumination
- Multi-connection turret allows attachment of most Storz, Olympus, Wolf and CMI fiber optic cables

- #49820 MFI Xenon Headlight System
- #49850 Xenon Mobile Stand



SFI Headlight Systems

The most advanced optics and material design culminate in an extremely lightweight, high performance headlight system. Exclusive condensing optics function to increase light intensity as spot size is reduced.

CL100 SFI Headlight System

- 100-Watt CL100 Surgical Illuminator
- 1000 hour average lamp life
- Xenon lamp provides high color temperature (6000K)
- System includes mobile stand

- #90137 CL100 Headlight System, 8' SFI fiber
- #90155 CL100 Headlight System, 12' SFI fiber

CL300 SFI Headlight System

- 300-Watt CL300 Illuminator
- 1000 hour average lamp life
- HID-MH lamp provides enhanced red spectral content
- System includes mobile stand

- #90130 CL300 SFI Headlight System, 8' SFI fiber
- #90157 CL300 SFI Headlight System, 12' SFI fiber

CL300 High Output Headlight System

Combining the durability of our MFI technology and the performance and comfort of our SFI technology

- Extreme performance
- CL300 Surgical Illuminator provides over 1000 hour average lamp life
- Standard FleeceComfort™ headband offers superior comfort
- Unique fused 2.8 mm shielded MFI fiber virtually eliminates all packing loss for maximum transmission efficiency
- System includes mobile stand

- #90140 CL300 High Output Headlight System



Optical Colposcope

- High-quality Leisegang™ optics with hand-polished glass and antiglare coating provide excellent depth of vision and clarity of focus
- Three-step magnification settings of 3.5x, 7.5x and 15x
- Simple optical measuring reticules
- Fine focus and vertical height adjustment

#88200 with Tilt Stand

#88460 with Swing Arm and Roller Base



Video Colposcope

Offer more complete care to your female patients with true video colposcopy that combines light, power and video for high-resolution imaging, superior diagnostics, and information management.


- The only direct-to-screen colposcope on the market gives you full-screen, high-resolution images
- Welch Allyn Solarc™ lighting generates 50% whiter, brighter light for true tissue color
- Store, retrieve, compare, annotate and transmit images with advanced image management
- Push button on-screen magnifications

#88000A with Vertical Arm Stand

#88001A with Vertical Arm Stand and Monitor

#89001A with Swing Arm Stand and Monitor

#88701 Image Management System for Colposcopy with Inkjet Printer

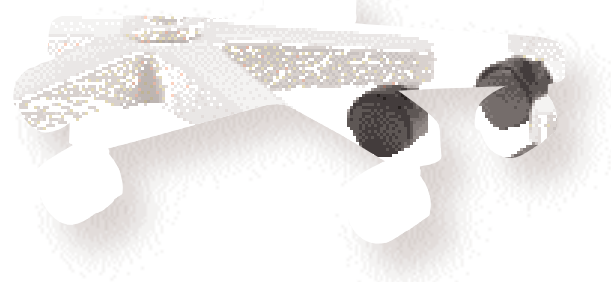

Exam Light III™


Flexible fiber optic lighting that lets you take brighter, whiter light anywhere without the heat.

- Durable fiber optic light pipe bends to almost any position and remains cool while transmitting illumination
- Small head leaves clear view of task areas
- Ideal for gynecological examinations and minor procedures
- Mounts easily on a wall, stand, or examination table
- Sealed focusing sleeve keeps optics free of dust
- Provides homogeneous spot and variable size

#48700 Exam Light III™

#48850 Mobile Stand

#48859 Universal Table Mount

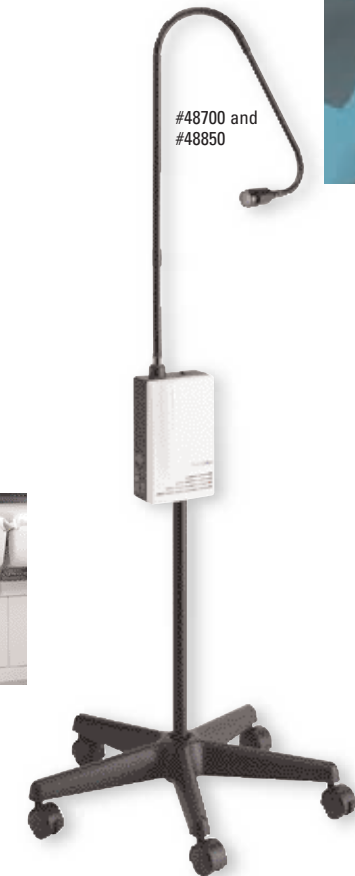

#58001

#78810

KleenSpec® Disposable Vaginal Specula and Illumination System

Conduct vaginal exams more completely and confidently with a vaginal specula that transmits brighter, whiter halogen lighting to the speculum tip and directly onto the cervix.

- Molded-in light carrier projects halogen lighting to the speculum tip and directly onto the cervix for clearer viewing and diagnosis
- A complete, permanent illumination system, including 6v halogen light handle, on/off switch, and transformer
- One-handed angulation and elevation provides multiple settings with positive-locking positions for easy operation
- Smooth, molded acrylic for maximum patient comfort
- Disposable for safety, economy and convenience

#78810 KleenSpec® Illumination System with Off/On Switch

#58001 KleenSpec® Medium Disposable Specula (Case of 100)

Electrosurgical System

Flexibility and efficiency in a mobile design.

- Generator and smoke evacuator in a single work station fits in even the smallest rooms
- Isolated power output and front LED allow precise power selection and delivery
- Microprocessor-controlled self-test mechanism, audible power alarm, and other safety features ensure accurate system operation
- Simple touch pad controls make it easy to operate and clean
- Full line of superior quality tungsten wire electrodes and filters available

#88530 Electrosurgical System with Cart and Smoke Evacuator


#88531 Electrosurgical Bundled Set with Cart

#88540 Electrosurgical Generator Only

#88550 Smoke Evacuator (Stand-alone)



#88530

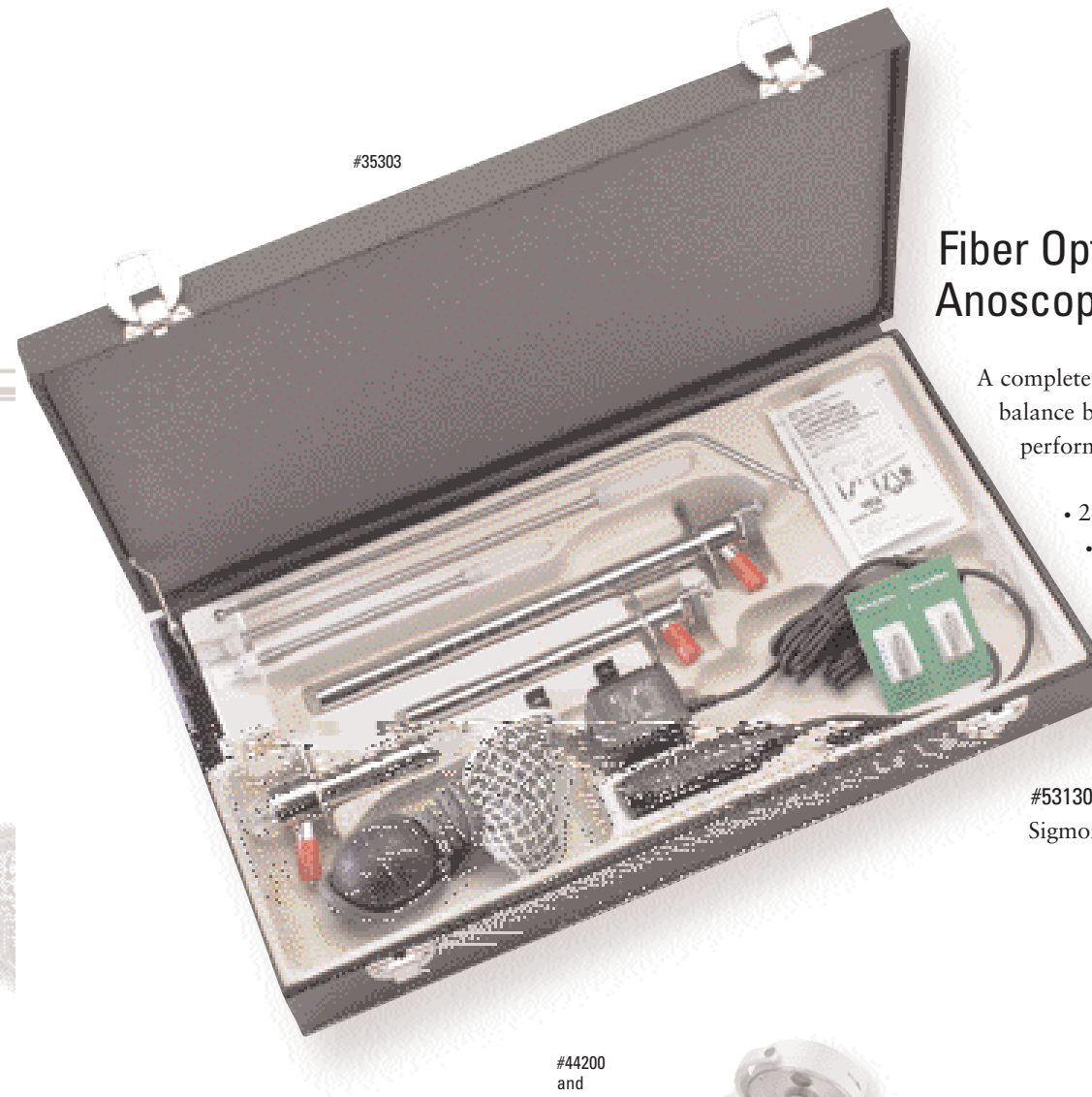

#35303

Fiber Optic Anoscope/Sigmoidoscope Set

A complete system that allows you to achieve the perfect balance between patient comfort and optimal clinical performance.

- 23mm x 7cm anoscope
- 15mm x 15cm pediatric, and 19mm x 25cm standard sigmoidoscopes
- Includes transformer, light handle with cord, complete insufflation bulb, two extra halogen lamps, suction tubes, and carrying case

#35303 Complete Fiber Optic Anoscope/Sigmoidoscope Set

#53130 Box of 25 KleenSpec® Disposable Sigmoidoscopes with Obturators

Cryosurgical System

All the features and components for expert care.

- Includes control console, five-foot probe line assembly, yoke, integrated probe storage area and 20-foot scavenging hose
- Lightweight design with convenient trigger lock and quick defrost for ease of operation and less finger fatigue
- Internal thermocouple allows precise monitoring of probe temperature for increased confidence
- Insulated probe barrel prevents freezing of adjacent tissue
- All probes and shaft assembly are sterilizable for patient safety

#88372 Cryosurgical System

#88373 Set of Three Cryosurgical Probes



#88372



#44200
and
#44250

LS-200™ Procedure Light


Bring light to all your procedures with high-quality, low-cost illumination.

- Superior maneuverability covers a wide range of positions without drifting or repositioning
- 24" working distance with 2,500 foot candles light output and 3,200K color temperature enhances visualization of the operative site and true tissue color
- Three-lamp design ensures continuous, shadow-free illumination
- 7" spot size with feathered edge definition and 12" depth of field reduces eye strain

#44200 LS-200™ Procedure Light

#44250 Mobile Stand

#44210 Wall Mount Bracket

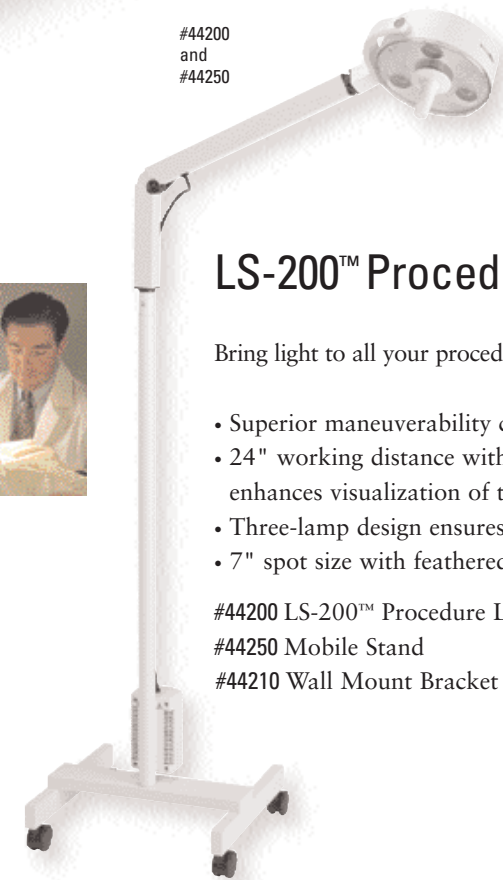

Consumables

Product Name	Part No.	Ref. Page	Product Affiliation
Eartips	#52434	14	524 Series KleenSpec® Disposable Otoscope Specula, 4 mm, 1,000/bag, 10 bags/case
	#52432	14	524 Series KleenSpec® Disposable Otoscope Specula, 2.5 mm, 1,000/bag, 10 bags/case
Thermometry Probe Covers	#05031-110	13	SureTemp® Disposable Probe Covers (10,000 covers)
	#05031-105	13	SureTemp® Disposable Probe Covers (5,000 covers)
	#05074-005	13	Braun Disposable Probe Covers for Pro3000 (case of 5,000)
Paper Supplies	#7052-25	10	Vital Signs Monitor Printer Paper, 1 case (6 rolls)
	#6200-40	5	Atlas™ Monitor Printer Paper, 1 case (25 rolls)
	#94016-0000	17	Paper for CP 10™ and CP 20™ 8.25" x 11" (1,000 sheets)
	#008-0040-98	6	Propaq® Monitor Standard Thermal Paper (10 rolls)
	#008-0178-00	6	Propaq® Monitor Fade Resistant Paper (10 rolls)
Batteries	#71960	14	Replacement Lithium Ion Battery
	#72200	14	Replacement Ni-Cad Rechargeable Battery (black) for Desk/Well Charger Handle
Lamps	#03100-U	14	Otosopes—20000, 25020, 21700, 20200, 3.5v Handle Adapter—73500 and Bivalve—26030
	#03800-U	14	PanOptic™ Ophthalmoscopes—11810, 11820
	#04200-U	23	Halogen Lamp for Exam Light II™—48600 and Exam Light III™—48700
	#04430-U	17	35w Halogen Lamp for LS-135 Exam Light—44300 and 44310
	#04450-U	17	50w Lamp for LS-150 Exam Light—44500 and 44510
	#04900-U	14	Ophthalmoscopes—11720, 11730, 11735
	#06000-U	19	Halogen HPX™ Lamp for Fiber Optic Laryngoscope Handles
	#06400-U	25	20w Lamp for LS-200™ Procedure Light—44200
	#07800-U	23-25	6v Halogen Lamp for Vag Spec Illuminator—78000, Sigmoidoscope, Anoscope Light Handles
	#09500-U	19-20	Solarc™ Lamp—49500, 46100
Disposable Sigmoidoscopes	#53130	(not pictured)	Individually Wrapped Disposable Sigmoidoscope Specula (19mmx25cm) with Obturators (25/box)
Disposable Anoscopes	#53110	(not pictured)	Individually Wrapped Disposable Anoscope Specula (19mmx10cm) with Obturators (25/box)
Electrodes	#45008-0000	17	ECG Electrodes, Resting for CP 10™ and CP 20™ (Box of 1,000)
	#6200-11	5	ECG Electrodes for Atlas Monitors (pack of 30)
Vaginal Specula	#58000	23	KleenSpec® 580 Series Disposable Vaginal Specula, Small (case of 100)
	#58001	23	KleenSpec® 580 Series Disposable Vaginal Specula, Medium (case of 100)
	#58004	23	KleenSpec® 580 Series Disposable Vaginal Specula, Large (case of 80)
	#58002	23	KleenSpec® 580 Series Disposable Barrier Vaginal Specula, Small (case of 100)
	#58003	23	KleenSpec® 580 Series Disposable Barrier Vaginal Specula, Medium (case of 100)
Miscellaneous	#11870	14	PanOptic™ Patient Eye Cups (package of 5)
	#52640	23	Disposable Sheaths for Exam Light II™ and Exam Light III™ (48600/48700)
	#58010	23	Disposable Sheaths for Vaginal Specula Illumination System (25 sheaths)

DURABLE ONE-PIECE CUFFS AND INFLATION SYSTEMS

MONITOR STYLE, ONE TUBE, SCREW CONNECTOR

Size	Circum.	Single	5 Pack*
Newborn	7.7-10.5 cm	5082-201-3	5082-81-3
Infant	9.8-13.3 cm	5082-202-3	5082-82-3
Small Child	12.4-16.8 cm	5082-203-3	5082-83-3
Child	15.8-21.3 cm	5082-204-3	5082-84-3
Child Print	15.8-21.3 cm	5082-204P-3	5082-84P-3
Small Adult	20.0-27.0 cm	5082-205-3	5082-85-3
Adult	25.3-34.3 cm	5082-206-3	5082-86-3
Adult Long	25.3-34.3 cm	5082-206L-3	5082-86L-3
Large Adult	32.1-43.4 cm	5082-207-3	5082-87-3
Large Adult Long	32.1-43.4 cm	5082-207L-3	5082-87L-3
Thigh	40.7-55.0 cm	5082-208-3	5082-88-3

Special Notes: To order multi-packs with one of each size cuff (excluding Child Print and Long Cuffs), order 5082-241-3.
 *To order cuffs with bayonet connectors (#5082-184), add "HP" suffix to part number (e.g., 5082-81-3HP).
 *To order cuffs with locking connectors (#5082-182), add "MQ" suffix to part number (e.g., 5082-81-3MQ).

MONITOR STYLE, TWO TUBES, SCREW CONNECTORS

Size	Circum.	Single	5 Pack*
Infant	9.8-13.3 cm	5082-202-4	5082-82-4
Small Child	12.4-16.8 cm	5082-203-4	5082-83-4
Child	15.8-21.3 cm	5082-204-4	5082-84-4
Child Print	15.8-21.3 cm	5082-204P-4	5082-84P-4
Small Adult	20.0-27.0 cm	5082-205-4	5082-85-4
Adult	25.3-34.3 cm	5082-206-4	5082-86-4
Adult Long	25.3-34.3 cm	5082-206L-4	5082-86L-4
Large Adult	32.1-43.4 cm	5082-207-4	5082-87-4
Large Adult Long	32.1-43.4 cm	5082-207L-4	5082-87L-4
Thigh	40.7-55.0 cm	5082-208-4	5082-88-4

Special Notes: To order multi-packs with one of each size cuff (excluding Child Print and Long Cuffs), order 5082-241-4.
 *To order cuffs with bayonet connectors (#5082-184), add "HP" suffix to part number (e.g., 5082-82-4HP).
 *To order cuffs with locking connectors (#5082-182), add "MQ" suffix to part number (e.g., 5082-82-4MQ).

SPHYG STYLE, ONE TUBE, LUER LOCK CONNECTOR

Size	Circum.	Single	5 Pack
Newborn	7.7-10.5 cm	5082-201-1	5082-81-1
Infant	9.8-13.3 cm	5082-202-1	5082-82-1
Small Child	12.4-16.8 cm	5082-203-1	5082-83-1
Child	15.8-21.3 cm	5082-204-1	5082-84-1
Child Print	15.8-21.3 cm	5082-204P-1	5082-84P-1
Small Adult	20.0-27.0 cm	5082-205-1	5082-85-1
Adult	25.3-34.3 cm	5082-206-1	5082-86-1
Adult Long	25.3-34.3 cm	5082-206L-1	5082-86L-1
Large Adult	32.1-43.4 cm	5082-207-1	5082-87-1
Large Adult Long	32.1-43.4 cm	5082-207L-1	5082-87L-1
Thigh	40.7-55.0 cm	5082-208-1	5082-88-1

Special Note: To order multi-packs with one of each size cuff (excluding Child Print and Long Cuffs), order 5082-241-1.

SPHYG STYLE, TWO TUBES, LUER SLIP CONNECTOR

Size	Circum.	Single	5 Pack
Infant	9.8-13.3 cm	5082-202-2	5082-82-2
Small Child	12.4-16.8 cm	5082-203-2	5082-83-2
Child	15.8-21.3 cm	5082-204-2	5082-84-2
Child Print	15.8-21.3 cm	5082-204P-2	5082-84P-2
Small Adult	20.0-27.0 cm	5082-205-2	5082-85-2
Adult	25.3-34.3 cm	5082-206-2	5082-86-2
Adult Long	25.3-34.3 cm	5082-206L-2	5082-86L-2
Large Adult	32.1-43.4 cm	5082-207-2	5082-87-2
Large Adult Long	32.1-43.4 cm	5082-207L-2	5082-87L-2
Thigh	40.7-55.0 cm	5082-208-2	5082-88-2

Special Note: To order multi-packs with one of each size cuff (excluding Child Print and Long Cuffs), order 5082-241-2.

INFLATION SYSTEMS, TWO TUBES WITH BULB AND VALVE

Size	Circum.	Single
Infant	9.8-13.3 cm	5082-222-1
Small Child	12.4-16.8 cm	5082-223-1
Child	15.8-21.3 cm	5082-224-1
Child Print	15.8-21.3 cm	5082-224P-1
Small Adult	20.0-27.0 cm	5082-225-1
Adult	25.3-34.3 cm	5082-226-1
Adult Long	25.3-34.3 cm	5082-226L-1
Large Adult	32.1-43.4 cm	5082-227-1
Large Adult Long	32.1-43.4 cm	5082-227L-1
Thigh	40.7-55.0 cm	5082-228-1

CONNECTORS AND ADAPTERS FOR DURABLE AND DISPOSABLE CUFFS

5082-161	Metal screw-type adapter with barbed end for Welch Allyn and Critikon monitors
5082-164	Plastic male screw-type connector with female luer slip for Welch Allyn and Critikon monitor cuffs
5082-165	Plastic male locking luer slip connector with barbed end for Datascope monitors
5082-168	Plastic female locking luer slip connector with barbed end for sphygmomanometers and Datascope monitor cuffs
5082-169	Metal female luer slip adapter with barbed end
5082-174	Plastic screw-type adapter with barbed end for Welch Allyn and Critikon monitors
5082-176	Plastic male luer slip connector with barbed end for neonatal cuffs
5082-178	Plastic male locking adapter for Welch Allyn Protocol and GE Marquette monitors
5082-182	Plastic female locking connector for Welch Allyn Protocol and GE Marquette monitor cuffs
5082-184	Plastic male bayonet connector for HP/Siemens, Spacelabs, and Datascope monitors
5082-185	Metal female bayonet adapter for HP/Siemens, Spacelabs, and Datascope monitor cuffs

DISPOSABLE ONE-PIECE CUFFS AND INFLATION SYSTEMS

MONITOR STYLE, ONE TUBE, SCREW CONNECTOR

Size	Circum.	Soft 5 Pack*	Vinyl 5 Pack*
Newborn	7.7-10.5 cm	5082-92-3	5082-91-8
Infant	9.8-13.3 cm	5082-92-3	5082-92-8
Small Child	12.4-16.8 cm	5082-93-3	5082-93-8
Child	15.8-21.3 cm	5082-94-3	5082-94-8
Small Adult	20.0-27.0 cm	5082-95-3	5082-95-8
Adult	25.3-34.3 cm	5082-96-3	5082-96-8
Adult Long	25.3-34.3 cm	5082-96L-3	5082-96L-8
Large Adult	32.1-43.4 cm	5082-97-3	5082-97-8
Large Adult Long	32.1-43.4 cm	5082-97L-3	5082-97L-8
Thigh	40.7-55.0 cm	5082-98-3	5082-98-8

Special Notes: To order multi-packs with one of each size cuff (excluding Long Cuffs), order 5082-241-7 for soft or 5082-241-19 for vinyl.
 *To order cuffs with bayonet connectors (#5082-184), add "HP" suffix to part number (e.g., 5082-91-3HP).
 *To order cuffs with locking connectors (#5082-182), add "MQ" suffix to part number (e.g., 5082-91-3MQ).

MONITOR STYLE, TWO TUBES, SCREW CONNECTORS

Size	Circum.	Soft 5 Pack*	Vinyl 5 Pack*
Infant	9.8-13.3 cm	5082-92-4	5082-92-9
Small Child	12.4-16.8 cm	5082-93-4	5082-93-9
Child	15.8-21.3 cm	5082-94-4	5082-94-9
Small Adult	20.0-27.0 cm	5082-95-4	5082-95-9
Adult	25.3-34.3 cm	5082-96-4	5082-96-9
Adult Long	25.3-34.3 cm	5082-96L-4	5082-96L-9
Large Adult	32.1-43.4 cm	5082-97-4	5082-97-9
Large Adult Long	32.1-43.4 cm	5082-97L-4	5082-97L-9
Thigh	40.7-55.0 cm	5082-98-4	5082-98-9

Special Notes: To order multi-packs with one of each size cuff (excluding Long Cuffs), order 5082-241-8 for soft or 5082-241-20 for vinyl.
 *To order cuffs with bayonet connectors (#5082-184), add "HP" suffix to part number (e.g., 5082-92-4HP).
 *To order cuffs with locking connectors (#5082-182), add "MQ" suffix to part number (e.g., 5082-92-4MQ).

MONITOR STYLE, ONE TUBE, LUER SLIP CONNECTOR

Size	Circum.	Soft 10 Pack
Neonate #1	3.3-5.6 cm	5082-101-1
Neonate #2	4.2-7.1 cm	5082-102-1
Neonate #3	5.4-9.1 cm	5082-103-1
Neonate #4	6.9-11.7 cm	5082-104-1
Neonate #5	8.9-15.0 cm	5082-105-1

Special Note: To order multi-packs with one of each size cuff, order 5082-241-9.

MONITOR STYLE, TWO TUBES, LUER SLIP CONNECTORS

Size	Circum.	Soft 10 Pack
Neonate #1	3.3-5.6 cm	5082-101-2
Neonate #2	4.2-7.1 cm	5082-102-2
Neonate #3	5.4-9.1 cm	5082-103-2
Neonate #4	6.9-11.7 cm	5082-104-2
Neonate #5	8.9-15.0 cm	5082-105-2

Special Note: To order multi-packs with one of each size cuff, order 5082-241-10.

SPHYG STYLE, ONE TUBE, LUER SLIP CONNECTOR

Size	Circum.	Soft 5 Pack	Vinyl 5 Pack
Newborn	7.7-10.5 cm	5082-91-1	5082-91-6
Infant	9.8-13.3 cm	5082-92-1	5082-92-6
Small Child	12.4-16.8 cm	5082-93-1	5082-93-6
Child	15.8-21.3 cm	5082-94-1	5082-94-6
Small Adult	20.0-27.0 cm	5082-95-1	5082-95-6
Adult	25.3-34.3 cm	5082-96-1	5082-96-6
Adult Long	25.3-34.3 cm	5082-96L-1	5082-96L-6
Large Adult	32.1-43.4 cm	5082-97-1	5082-97-6
Large Adult Long	32.1-43.4 cm	5082-97L-1	5082-97L-6
Thigh	40.7-55.0 cm	5082-98-1	5082-98-6

Special Note: To order multi-packs with one of each size cuff (excluding Long Cuffs), order 5082-241-5 for soft or 5082-241-17 for vinyl.

SPHYG STYLE, TWO TUBES, LUER SLIP CONNECTOR

Size	Circum.	Soft 5 Pack	Vinyl 5 Pack
Infant	9.8-13.3 cm	5082-92-2	5082-92-7
Small Child	12.4-16.8 cm	5082-93-2	5082-93-7
Child	15.8-21.3 cm	5082-94-2	5082-94-7
Small Adult	20.0-27.0 cm	5082-95-2	5082-95-7
Adult	25.3-34.3 cm	5082-96-2	5082-96-7
Adult Long	25.3-34.3 cm	5082-96L-2	5082-96L-7
Large Adult	32.1-43.4 cm	5082-97-2	5082-97-7
Large Adult Long	32.1-43.4 cm	5082-97L-2	5082-97L-7
Thigh	40.7-55.0 cm	5082-98-2	5082-98-7

Special Note: To order multi-packs with one of each size cuff (excluding Long Cuffs), order 5082-241-6 for soft or 5082-241-18 for vinyl.

INFLATION SYSTEMS, TWO TUBES WITH BULB AND VALVE

Size	Circum.	Soft Single	Vinyl Single
Infant	9.8-13.3 cm	5082-232-1	5082-232-6
Small Child	12.4-16.8 cm	5082-233-1	5082-233-6
Child	15.8-21.3 cm	5082-234-1	5082-234-6
Small Adult	20.0-27.0 cm	5082-235-1	5082-235-6
Adult	25.3-34.3 cm	5082-236-1	5082-236-6
Adult Long	25.3-34.3 cm	5082-236L-1	5082-236L-6
Large Adult	32.1-43.4 cm	5082-237-1	5082-237-6
Large Adult Long	32.1-43.4 cm	5082-237L-1	5082-237L-6
Thigh	40.7-55.0 cm	5082-238-1	5082-238-6